


# Mold Composite Auxiliaries Directly on Teeth

## EASY, EFFECTIVE & ECONOMICAL

- Molds are translucent and flexible to facilitate composite loading, curing, and disengagement.
- A mold is filled with a light cure composite (Go To, TKO, or Flowtaint recommended) then delivered to the tooth and cured.
- The desired auxiliary is delivered, formed, and functional immediately.
- Value priced to be single-use disposable.

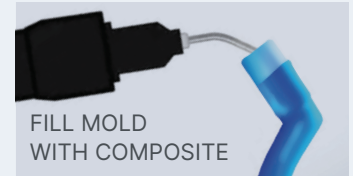

MOLDS MAY BE DELIVERED WITH BONDER HANDLE (#BHA) OR ANY HEMOSTAT.

## DIFFERENT MOLDS FOR DIFFERENT APPLICATIONS

ONLY FROM  **RELIANCE®**



### HABIT CONTROL

Inhibits thumbsucking and tongue thrust tendencies.


Habit control #RMH



### CUSPID RAMP

Bite opening ramp applied to the cuspid—5.0 mm.


Cuspid ramp #RMC

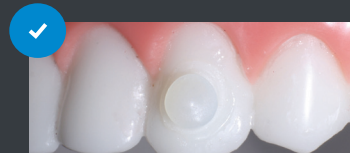

### BUTTON

Low profile and smooth. Facilitates elastomeric or ligature engagement.


Button #RMBU

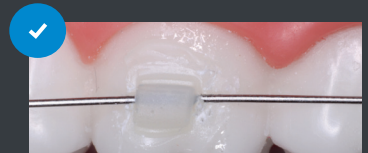

### MOLAR TUBE

Molds a temporary tube. Accepts wire up to 020.

Molar Tube #RMT

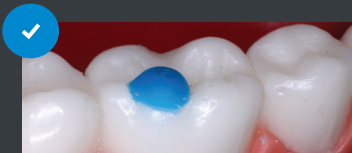

### MOLAR PAD

Open the bite with posterior occlusal buildups. TKO or Ultra Band Lok recommended.


Molar Pad #RMMP



### BITE RAMP

Bite opening turbo formed on the lingual of upper incisors. 2 sizes—3.0 mm or 5.0 mm.

Bite ramp 3mm #RMBR3  
Bite ramp 5mm #RMBR5

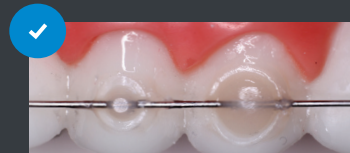

### WIRE HOLDER

Bonds lingual retainers or other wires consistently with a smooth and low profile finish. 2 sizes—3.0 mm or 4.5 mm.

Wire bonder large #RMWBL  
Wire bonder small #RMWBS

### MOLD KIT

(2) Bonder Handles  
(108) Molds; 12 Per Type

Mold Kit #RMK

### INDIVIDUAL MOLDS

(25) Molds Per Package

Bonder Handle Sold Separately  
#BHA

 RELIAMOLD INSTRUCTIONAL VIDEO AVAILABLE AT [RELIANCEU.TECH](http://RELIANCEU.TECH)