


* Ortho KIT

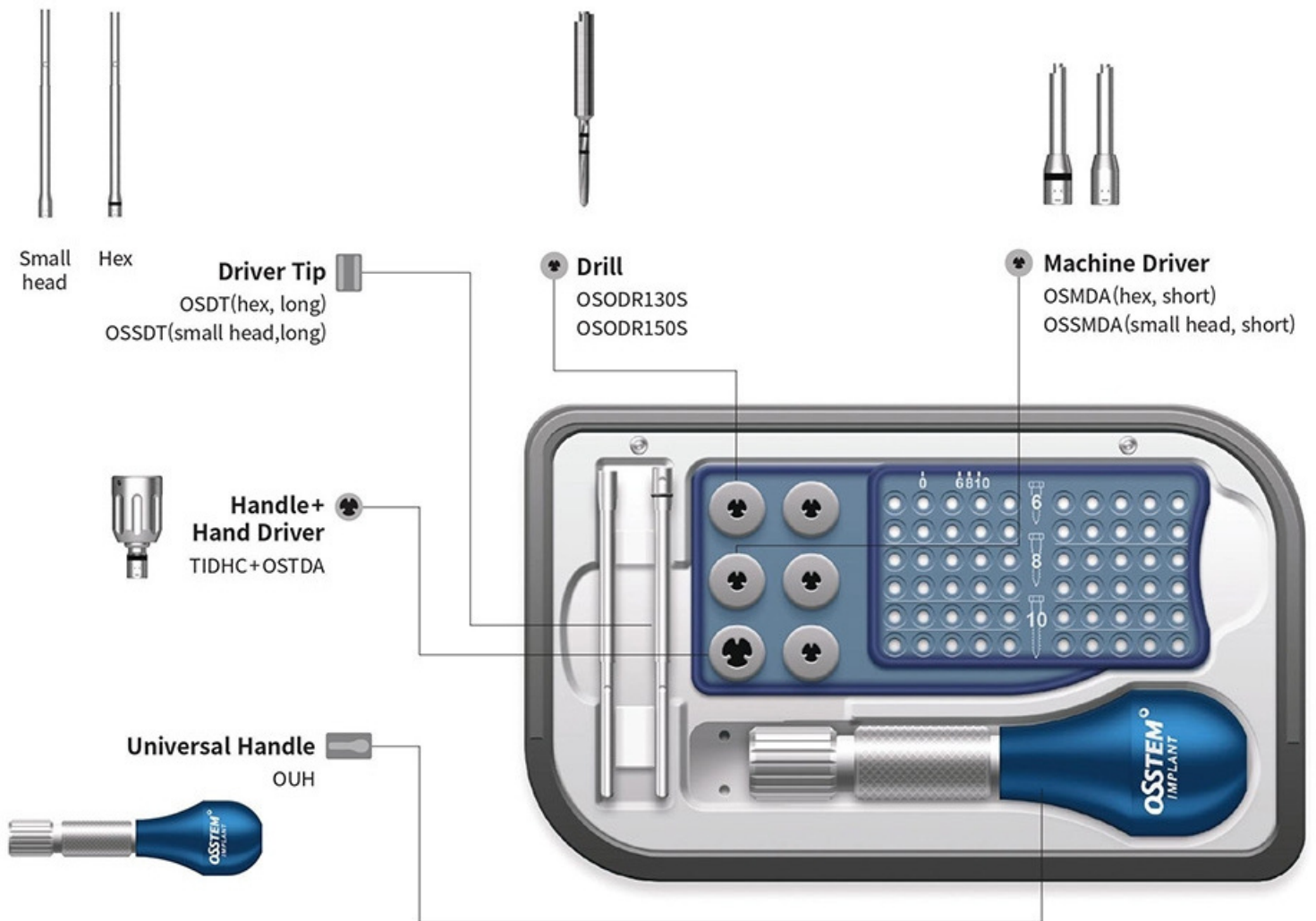


Order NO.

OOKS

* Drill

- Use drill by connecting with handpiece. (engine)
- Ø1.0 drill: use before Ø1.2/1.4 screw insertion.
- Ø1.3 drill: use before Ø1.6 screw insertion.
- Ø1.5 drill: use before Ø1.8 screw insertion.
- Recommend speed of drill: 800 rpm. (high speed)
- Before screw insertion recommend to the eliminated cortical bone only. (if the cortical bone were very thick, drilling equal length as screws)
- Ø1.0 drill is not included in the Kit.

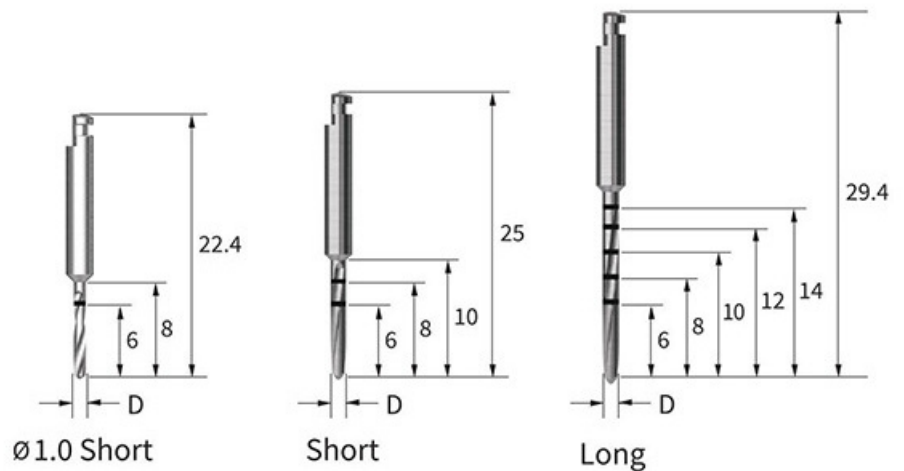


L \ D	Ø1.0	Ø1.3	Ø1.5
Short	OSODR100S	OSODR130S	OSODR150S
Long	-	OSODR130C	OSODR150C



* Universal Handle

- The universal handle can be used by connecting with the driver tip to screw insertion.
- The scribed grip to support manual insertion.

Order NO.	OUH
-----------	-----

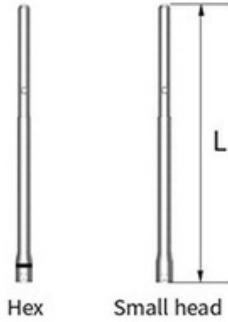

* Driver Tip

- Driver tip use to screw insertion by connecting with Universal handle.
- Hex driver designed to connect with normal head screw.
- Small head driver designed to connect with small head screw.

L Type	Hex	Small Head
Short(L)	OSDTS(45)	OSSDTS(45)
Long(L)	OSDT(67)	OSSDT(67)



* Driver Handle

- Driver handle use to screw insertion by connecting with a Hand driver.


Order NO.	TIDHC
-----------	-------



* Hand Driver

- Hand driver designed to screw insertion by connecting with Driver handle or ratchet wrench.
- Two types of Hand driver to use different size of the screw head.
- Hex driver for normal size of the screw head, and small driver for the small size of the screw head.
- Small head hand driver is not included in the Kit.

Hex	Small Head
OSTDA	OSSTDA

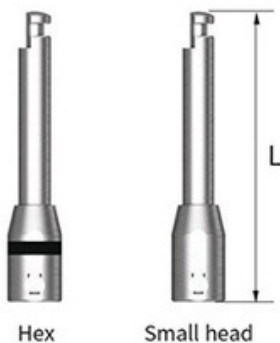

* Machine Driver

- Machine driver designed to screw insertion by connecting with engine.
- Two types of hand driver to use different size of the screw head.
- Hex Machine driver for normal size of the screw head.
- Small head Machine driver for the small size of the screw head.

L Type	Hex	Small Head
Short(L)	OSMDA(21.4)	OSSMDA(21.4)
Long(L)	OSMDB(31.4)	OSSMDB(31.4)

*** MINI SCREW TOOL**



Surgical Tray [801-61]



Handle Driver [801-62]



Small Hand Driver [801-66]



Driver Shaft [801-63]
(All kinds of Handle Driver can be used)



Long Tip [801-64]



Short Tip [801-65]



Drill Tip [801-67]

Mini Screw Tool

Type	Order NO.
Mini Screw Tool Kit	800-00
Surgical Tray	801-61
Handle Driver	801-62
Driver Shaft	801-63
Long Tip	801-64
Short Tip	801-65
Small Hand Driver	801-66
Drill Tip	801-67